

ROUND, BICONVEX) & ANATOMICAL GLUTEAL IMPLANTS



SEBBIN PARIS
ESTHETIQUE & RECONSTRUCTION

Recreating harmonious bodies for a new vision of rebirth.

Anatomical gluteal implants




| Reference | Volume (mL) | Width (mm) | Height (mm) | Projection (mm) |
|-----------|-------------|------------|-------------|-----------------|
| LS 05 185 | 185 | 97 | 139 | 30 |
| LS 05 260 | 260 | 107 | 148 | 32 |
| LS 05 335 | 335 | 117 | 162 | 36 |
| LS 05 400 | 400 | 129 | 174 | 42 |

This shape is recommended for patients when:

- The majority of the gluteal volume is in the upper part of the buttock,
- The gluteal muscle is long,
- The buttocks form a high square,
- The patients are tall.

It mainly allows to increase the volume of the lower part of the buttock.

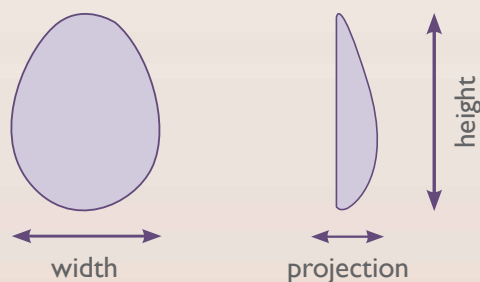


Round gluteal implants

| Reference | Volume (mL) | Base (mm) | Projection (mm) |
|-----------|-------------|-----------|-----------------|
| LS 04 200 | 200 | 108 | 30 |
| LS 04 235 | 235 | 113 | 33 |
| LS 04 270 | 270 | 120 | 34 |
| LS 04 300 | 300 | 125 | 36 |
| LS 04 330 | 330 | 127 | 38 |
| LS 04 360 | 360 | 128 | 42 |
| LS 04 400 | 400 | 132 | 43 |
| LS 04 440 | 440 | 136 | 44 |

This implant is recommended in the following indications:

- If the majority of the gluteal volume is in the lower part of the buttock,
- If the gluteal muscle is short,
- If the buttocks form a small square.

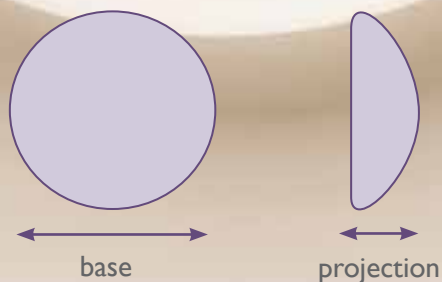


Biconvex gluteal implants



| Reference | Volume (mL) | Base (mm) | Projection (mm) |
|-----------|-------------|-----------|-----------------|
| LS 06 370 | 370 | 127 | 52 |
| LS 06 410 | 410 | 130 | 53 |
| LS 06 480 | 480 | 136 | 56 |
| LS 06 530 | 530 | 143 | 59 |

- The biconvex shape answers the same indications as the round shape but ensures more projection for an identical base.
- Moreover, in the event of a shift from deep face to superficial face, the biconvex implant being in all points symmetrical, makes this displacement indiscernible for the patient.

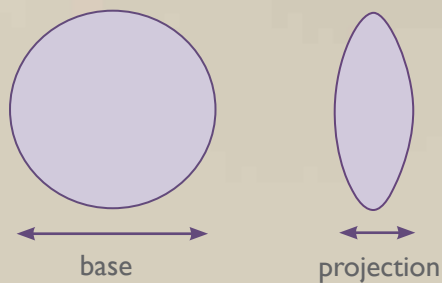

Features



Groupe SEBBIN is one of the first European manufacturers to propose this implantable medical device.

The gluteal implants prefilled with silicone gel are more flexible, easier to introduce and offer a better projection compared to the silicone elastomer monobloc implants.

During the introduction, the implant can be strongly stressed; thanks to their specific very high cohesive silicone gel, our gluteal implants return rapidly to their initial shape.

Our sterile gluteal implants, for single use, are sold by unit. They can be manufactured upon request with a textured envelope under the reference LST.

The dimensions of our round, biconvex and anatomical implants can differ slightly from the data mentioned in this brochure.



Our commitments



CONTINUITY

Production know-how mastered for almost 30 years.

RELIABILITY

Biological and mechanical evaluation of end products.
Visual controls on 100 % of implants at each step of the manufacturing process.

TRACKING

Each implant is identified by a unique serial number allowing tracking from the raw materials used to its implantation.

PEACE OF MIND

Highly qualified personnel working in clean room.

This document is dedicated to health care professionals. According to the Medical Device Directive 93/42/EEC, the gluteal implants are Class IIb medical devices manufactured by Groupe SEBBIN and are intended to be used in the context of a plastic, reconstructive and cosmetic surgery. These devices are CE marked by the notified body MDC number 0483. Please read carefully the instructions for use before using.

MANUFACTURER:



GROUPE SEBBIN SAS
39, Parc d'Activités des Quatre Vents
95650 Boissy l'Aillerie - France
Tel.: +33 1 34 42 13 28
Fax: +33 1 34 42 16 88
www.sebbin.com

INTERNATIONAL OPERATIONS:

Sebbin Benelux SPRL:
Tel.: +32 3 451 26 81 - Fax: +32 3 451 26 89
Email: info@laboratoires-sebbin.be
Sebbin Deutschland GmbH:
Tel.: +49 2102 70 88 00 0 - Fax: +49 2102 70 88 00 99
Email: service@sebbin-implantate.de
Groupe Sebbin SAS, France:
Tel.: +33 1 34 42 13 28 - Fax: +33 1 34 42 16 88
Email: serviceclients@laboratoires-sebbin.fr
Sebbin Iberica Productos Medicos, S.L.:
Tel.: +34 93 159 3003 - Fax: +34 93 676 3053
Email: info@sebbin.es
Sebbin Suisse Sarl:
Tel.: +41 43 336 72 30 - Fax: +41 43 336 72 36
Email: info@sebbin.ch
Sebbin UK Ltd:
Tel.: 0800 652 4028 - Fax: +4433 0088 9575
Email: custserv@sebbin.uk